

## LUNA™ 8-Channel Slides

L72001 (1 box of 50 slides)  
L72002 (10 boxes of 50 slides)  
L72003 (10 boxes of 50 slides)  
Sterile - gamma-irradiated)

## Product Description

Materials	PMMA
Dimensions (W x D x H)	75 mm x 25 mm x 2.1 mm
Number of chamber	8 chambers
Chamber volume	10 $\mu$ L

LUNA™ 8-Channel Slides are individually wrapped disposable slides used with the LUNA-FX7™ Automated Cell Counter.



## Storage

Room temperature

- Sample throughput: Up to 8 samples
- Sample loading volume: 10  $\mu$ L/chamber
- Maximum analysis volume: 0.5  $\mu$ L/chamber

## Directions for Use

1. Prepare sample to count with or without staining dye.
2. Tear off the top of the wrap and gently take out the slide.
3. Load 10  $\mu$ L of a sample into each chamber.
4. Multichannel dispensing can be used with a multichannel pipette.
5. Apply the slide to the LUNA-FX7™ Automatic Cell Counter.



## HEADQUARTERS

FL 3  
28 Simindaero 327beon-gil, Dongan-gu  
Anyang-si, Gyeonggi-do 14055  
South Korea

Tel: +82 (31) 478-4185

## USA

7700 Little River Turnpike STE 207  
Annandale, VA 22003  
USA

Tel: +1 (703) 622-4660, +1 (703) 942-8867

## EUROPE

1 allée Lavoisier 59650 Villeneuve d'Ascq  
France

Tel: +33 (0)3 74 09 44 35

[www.logosbio.com](http://www.logosbio.com)

LBSM-RD-PI-L8-001 Rev.1



## Product Care

- ! Do not press down a slide with or without wrapping.
- ! When handling, hold the edges of the wrapping, not the slide itself.
- ! Widen the opening and gently remove the slide to prevent possible scratches.
- ! Avoid fingerprints on the chamber area.
- ! Do not load any solvents buffer for counting.

## Disclaimer

This product is for research use only.  
Please consult the material safety data sheet for information regarding hazards and safe handling practices.

Additional information is available on our website [www.logosbio.com](http://www.logosbio.com).  
© Aligned Genetics, Inc. 2021. All rights reserved.

## LUNA™ 8-Channel Slides

L72001 (1 box of 50 slides)  
L72002 (10 boxes of 50 slides)  
L72003 (10 boxes of 50 slides)  
Sterile - gamma-irradiated)

## Product Description

Materials	PMMA
Dimensions (W x D x H)	75 mm x 25 mm x 2.1 mm
Number of chamber	8 chambers
Chamber volume	10 $\mu$ L

LUNA™ 8-Channel Slides are individually wrapped disposable slides used with the LUNA-FX7™ Automated Cell Counter.



## Storage

Room temperature

- Sample throughput: Up to 8 samples
- Sample loading volume: 10  $\mu$ L/chamber
- Maximum analysis volume: 0.5  $\mu$ L/chamber

## Directions for Use

6. Prepare sample to count with or without staining dye.
7. Tear off the top of the wrap and gently take out the slide.
8. Load 10  $\mu$ L of a sample into each chamber.
9. Multichannel dispensing can be used with a multichannel pipette.
10. Apply the slide to the LUNA-FX7™ Automatic Cell Counter.



## HEADQUARTERS

FL 3  
28 Simindaero 327beon-gil, Dongan-gu  
Anyang-si, Gyeonggi-do 14055  
South Korea

Tel: +82 (31) 478-4185

## USA

7700 Little River Turnpike STE 207  
Annandale, VA 22003  
USA

Tel: +1 (703) 622-4660, +1 (703) 942-8867


## EUROPE

1 allée Lavoisier 59650 Villeneuve d'Ascq  
France

Tel: +33 (0)3 74 09 44 35

[www.logosbio.com](http://www.logosbio.com)

LBSM-RD-PI-L8-001 Rev.1

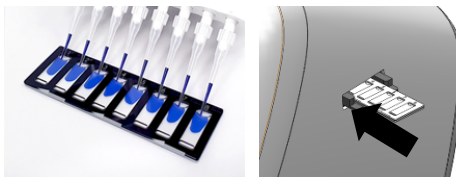

## Product Care

- ! Do not press down a slide with or without wrapping.
- ! When handling, hold the edges of the wrapping, not the slide itself.
- ! Widen the opening and gently remove the slide to prevent possible scratches.
- ! Avoid fingerprints on the chamber area.
- ! Do not load any solvents buffer for counting.

## Disclaimer

This product is for research use only.  
Please consult the material safety data sheet for information regarding hazards and safe handling practices.

Additional information is available on our website [www.logosbio.com](http://www.logosbio.com).  
© Aligned Genetics, Inc. 2021. All rights reserved.